

**Mounting.
Anything. Anywhere.**
...for more than 20 years

medically certified



Your competent mounting partner

Benefit from Rehadapt as your competent mounting partner for **medical devices** and more!

For two decades, Rehadapt has been globally recognized for mounting systems that securely **position assistive technologies** for optimal user experience.


Our extensive experience in the medical device sector makes us the ideal partner for **all types of professional mounting requirements**.

Rehadapt is a dedicated team at your side that will always find your very perfect mounting solution.

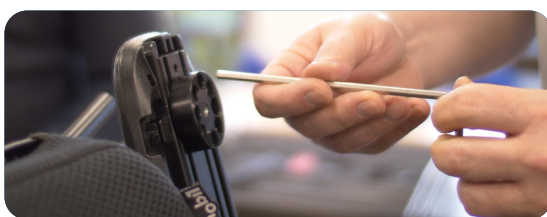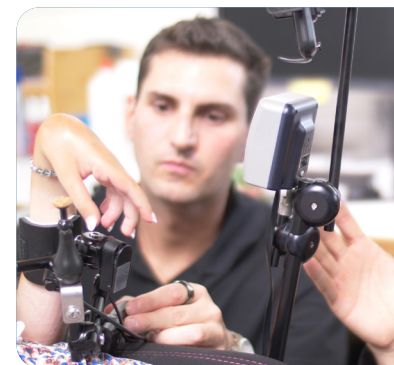

The individual needs of our users and regulatory requirements have driven the development of the most **comprehensive modular system** for mounting solutions worldwide.

Our components can be **combined almost playfully**, much like building blocks. It's the perfect choice for those who want to explore new and creative solutions, and it's always in compliance with medical device standards!




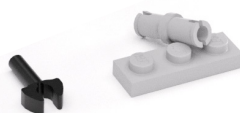


Das Rehadapt-System ist bestens geeignet für:

- ✓ Environmental controls
- ✓ Speech devices
- ✓ iPads
- ✓ Smartphones
- ✓ Monitors
- ✓ ventilation hose
- ✓ ventilation equipment
- ✓ stationary pulse oximeter
- ✓ blood pressure monitoring equipment
- ✓ long-term ECGs
- ✓ Input devices
- ✓ Storage systems
- ✓ Tool holders
- ✓ Lighting
- and much more...

In folgenden Bereichen:

- ✓ in medical sector
- ✓ in laboratory settings
- ✓ in the rehabilitation sector
- ✓ in the care environment
- ✓ in workshops
- ...and always MDR-compliant!

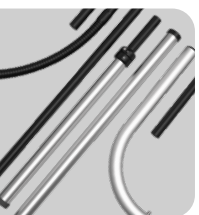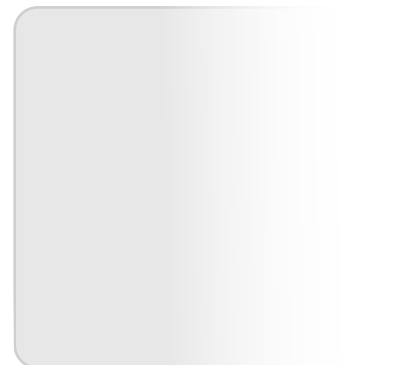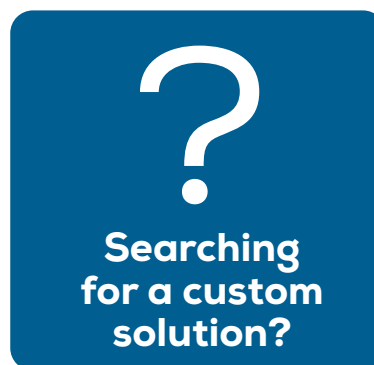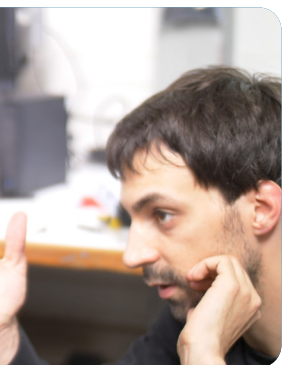

Would you like to learn more about Rehadapt mounting systems or **are you searching** for a custom clamp for your product? Then please feel free to contact us.


Our team is looking forward to hearing from you! We mount anything and anywhere - reliably and professionally for more than 20 years!




How to tackle complex challenges?

With a system!

With every new mounting challenge, essentially the same questions arise: Which device do I want to fasten? Where should it be positioned? Where is the nearest mounting point? And how do I bridge the distance from the mounting point to the desired position? Our **well-thought-out component system** leaves no limits to your creativity.

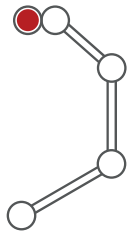

Device adapter

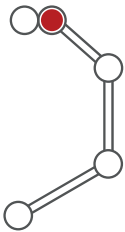


Quick-release socket



Joints



Tubes

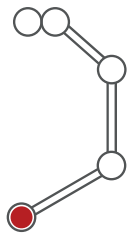

Mounting basis



Kits




Tools




Replacement parts

Our decades of experience in the development, manufacturing, and application of our components have led to a wide range of **configurable mounting solutions**. This encompasses all areas for the positioning of assistive technologies and medical devices.

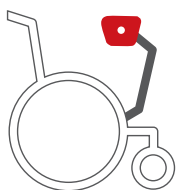

Wheelchair mounts

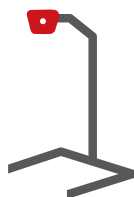


Floor stands




Table mounts




Wall mounts

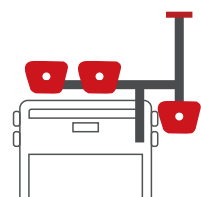

Out-patient care

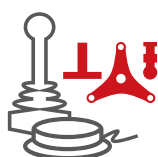

Control mounts

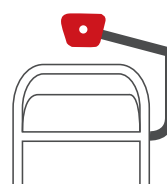


Bed mounts



Workshop mounts

Our „**Universal Device Socket „ (UDS)** includes our quick re-release mechanism which allows for the safe device attachment within seconds. The footprint on the device is minimal, yet great sturdiness is achieved. This allows any device to **easily switch between different applications** and their respective mounts using an appropriate adapter.

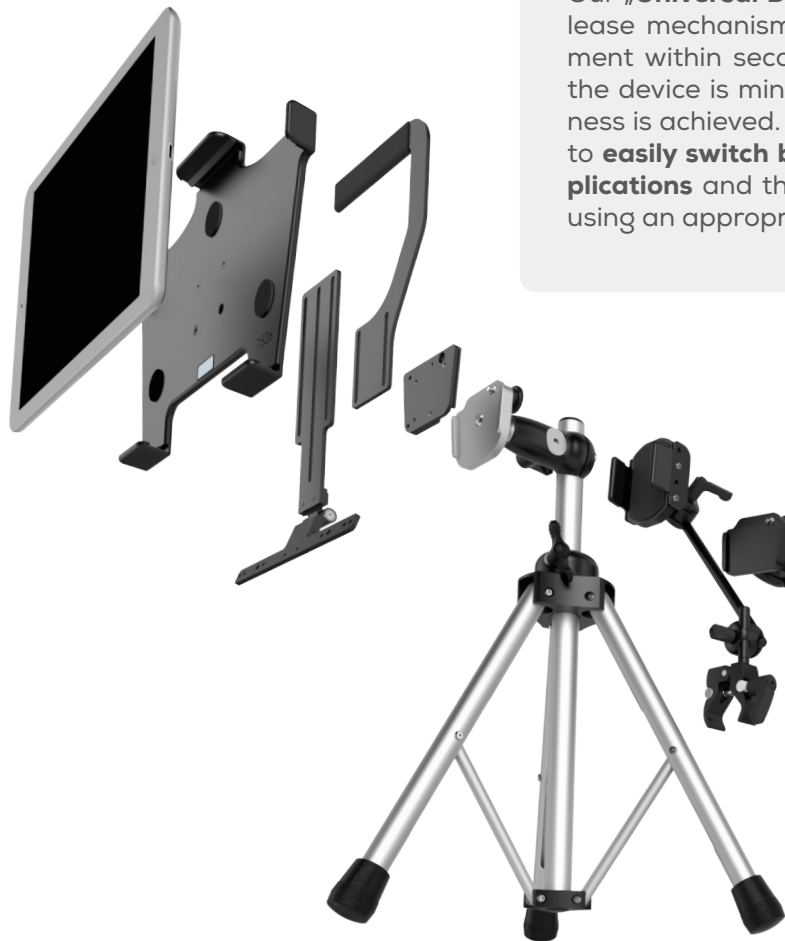

**360° of
mounting**

Do you think your user could benefit from additional mounts for their devices?


How about a speech generating device for example? Is there a need of a floorstand in their living room or a table stand at school? How do they communicate while travelling?

Open your mind for **a holistic view on your client's needs** and full mounting assessments.

There are **no limits to the extendability** of our mounting solutions. Whether on the floor, on a table, or on a frame profile, the Rehadapt system can be precisely adapted for any use case.




Mounting-Base

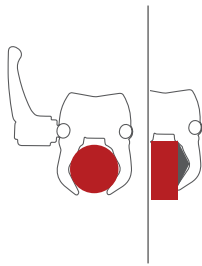

MiniClamp

Quick to mount without taking up much space. Featuring rubberized clamps and additional inserts for rectangular and round tubes. Ideal for lightweight devices such as switches.

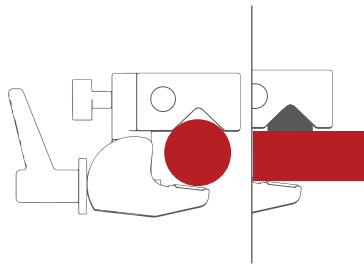


SuperClamp

Designed for a fast and highly stable mounting solution. Can be attached to both round and square elements, such as tables or armrests.




Wall mounts

Featuring large plates for optimal weight distribution. Can be mounted vertically or parallel to the wall. Available with a complete 3D arm or as a parking space for an existing mount.

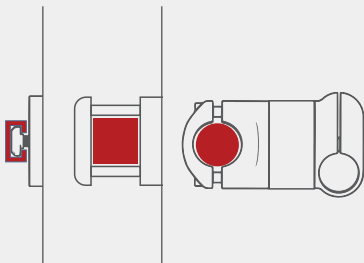


Frame clamps

Frame clamps are designed to precisely match the respective frame profile. Ideally suited for permanently installed mounting solutions.

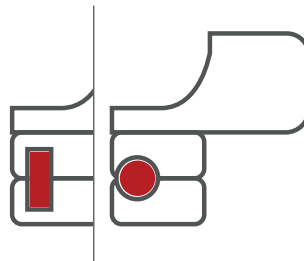

Bed clamps

Clamps directly onto the slatted frames of well-known nursing bed manufacturers. Moves the entire mounting solution accurately with the movement of the headboard.

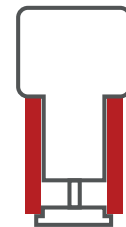

Insert sleeves

Insert sleeve for attaching various mounting solutions in the field of outpatient care. Installed directly on the patient lifter receptacle of nursing beds.

Each mounting base consists of at least two parts. There is always a frame clamp and a tube socket. In addition, there are extensions in the middle: distance spacers, angles, or a coupler.



Round



Rectangular



Rails



Holes



closed



semi-open



on flat surfaces



on tubes

Device adapters & sockets



Quickly hook them in and out of the universal device socket.




Insert them into the spigot link system.



Screw them onto the MagicArm.


Joints

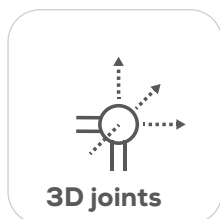

Full 360° adjustment in all dimensions



Swing your mounts between different positions

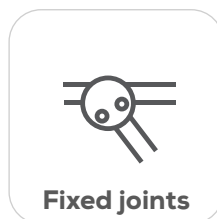

Connect tubes at a fixed angle



Move a device easily out of the way



Activation with a simple push button

Tubes



Diameters



“With Rehadapt you are not on your own.
Through our wide range of services, you’ve got a
dedicated team at your side.”



Obtain professional tools and expert training to serve your clients with profound positioning and mounting skills. We are happy to share our knowledge and transfer mounting expertise to the field through **on-demand videos, online trainings and workshops.**



Your communication device deserves the perfect mount and you deserve superior service. We analyze your **individual** needs to provide you with the best **mounting solution.** For free!



Rehadapt wants to **increase the lifetime of products** and provide support for its clients even after the 5 year limit which is specified by law through the Medical Device Regulation. Maintenance in the RSC facilitates the usage of products until the end of their true service time. Furthermore we offer structured **maintenance documents** for correctly reviewing Rehadapt mounting solutions.

Learn more on rehadapt.com and feel free to get in touch with us!



Rehadapt Engineering GmbH & Co. KG
Heinrich-Hertz-Str. 104
34123 Kassel, Germany
Fax +49 561 220 717 0
Fon +49 561 220 717 99
officeteam@rehadapt.com

Don't hesitate and give us a call. You can download this flyer or our current product- and service catalogs here:



www.rehadapt.com/downloads