



# REHADAPTER

The technological progress of consumer devices has made tablets like the iPad ready for serious AAC tasks. Let us show you how an iPad can be your next speech generating device.


Welcome ...to the Rehadapter.

### THE REHADAPTER IS READY FOR EVERY AAC TASK - WHENEVER YOU NEED IT

The Rehadapter seamlessly transforms the iPad into a dedicated AAC device. It has everything you would expect from the top performing AAC devices on the market: clear and loud sound, connections for eye trackers, keyguards, switches, a kickstand for quick positioning, maximum portability, and a built-in universal device connector for the most advanced mounting system on the market. And all this comes within one protective enclosure.

Combined with the iPad's bandwidth of apps, the Rehadapter offers endless possibilities for AAC users. It is the only case on the market to offer this level of flexibility, technology, and ease of use. The Rehadapter is designed for all-day communication while at home, at school or travelling.

Choose the Rehadapter, because now more than ever.... we hold your independence.

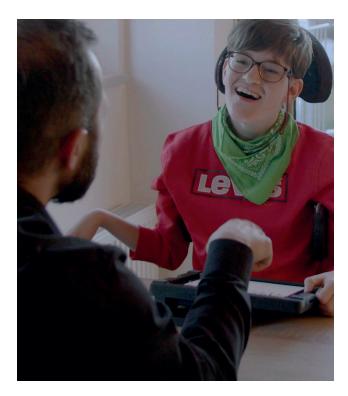



#### **ALL-IN-ONE**

The Rehadapter transforms your iPad into a dedicated AAC device by combining all key features in one single housing!













#### **FUTURE-PROOF & MODULAR**

By design the Rehadapter is future-proof when it comes to extendability. All access methods you might want to use in future are covered. It's the first case on the market with this broad connectivity. And what if you need to replace your iPad? We will easily and quickly be able to accomodate for the next generations of iPads or other selected tablets.

"We are proud that the Rehadapter turns an ordinary iPad into a full-fledged AAC tool for helping people with special needs finally hassle-free and more flexible then ever while being economical!"



**ULI EHLERT**CEO Rehadapt

# Your new AAC hub equipped with all features needed



The Rehadapter is a protective case designed for extensive daily usage. It has a sturdy design and is made of **reinforced plastic**, delivering ample protection. The iPad is surrounded by rubber bumpers that absorb impacts from all sides.



The **magnetic charging cable** ensures not only self-alignment when connecting, but also makes sure the Rehadapter cannot be pulled or dragged by its cable.



All in all it's a trusty companion which meets all the requirements for a **medical product** and is **scratch**, **shock and drop resistant**.



High-quality sound is achieved by one large integrated speaker located on the back of the unit. With its **voice-optimized** loud and clear sound the Rehadapter is classroom-ready.

The rubberized volume buttons are always accessible and the angle of the speaker has been optimized for **maximum performance** in all positions: whether the Rehadapter is flat on a table, upright, or held by a Rehadapt mounting system.





The **battery lasts** twice as long as the iPad. This means you can use it without having to worry about it throughout the day.



AAC solutions only work well when the client's access needs are perfectly matched. Never before has an adapted iPad offered so many options.

#### **UNIQUE KEYGUARD SYSTEM**

The Rehadapter features an innovative keyguard system unmatched in the market: the keyguards can be **clicked in place** without the need of any additional holder, velcro or moving piece.

This unique system can be set to **three different levels of pressure** to better fit the client's needs: easy toolless, semi-permanent, or permanent (requiring tools for removal).



The keyguards are made from 5 mm thick plexiglass and can be supplied off-the-shelf or customized. They are milled instead of lasered to make sure the edges are **smooth and comfortable** for the user.

#### TWO SWITCH ACCESS

The Rehadapter features **integrated 3.5 mm jacks** to connect up to two switches for screen scanning. No need for additional, expensive and power consuming Bluetooth switches that might leave you with an unreliable connection.



#### **EYETRACKER, HEAD MICE & MORE**

The Rehadapter features **two attachment points** for additional access equipment such as eye trackers or head mice devices. The attachment points - located at the top and bottom of the device - offer **unlimited flexibility** for attaching additional equipment.



The Rehadapter includes an **integrated kickstand** offering continuous adjustment to an inclination angle of up to 45°. When lying flat on the table the Rehadapter's wedge shape allows for an ergonomic 10° angle access. The kickstand and all contact points are rubberized. The kickstand is made of aircraft-grade aluminum and can be removed if not needed.

Rehadapt's universal device socket (UDS) is an established standard in the AAC field and the Rehadapter naturally includes an **adapter plate**. This makes it immediately ready to be connected with the most comprehensive mounting ecosystem on the market.





The Rehadapter features a sophisticated **shoul-der strap system** that will never become twisted. The Rehadapter can also lie flat on a surface without wobbling even when the strap is attached. The **handle** is wide and rubberized for comfortable grip. That makes it easy to carry the Rehadapter around or to stow it into a bag.

#### **SAFETY BOLT**


The case is closed fasteners and n in. This provides it

#### KEYGUARD

Due to the vast amount of AAC apps on the market keyguards can be customized for any particular need. Over time keyguards for all standard apps will be supplied off-the-shelf.

#### **FRONT FRAME**

The Rehadapter is future-proof for any iP upcoming to the market. The only part wh has to be newly engineered is the according front frame. So for any new iPad there will be fitting Rehadapter in no time.

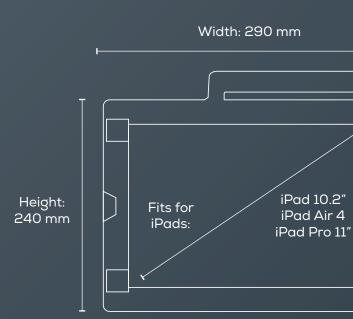

# REHADAPTER



d with mechanical ot simply clicked-maximum safety.

ad ich

ng e a

ב

## DEDICATED VOLUME CONTROL


The Rehadapter features custom volume buttons separated from the iPad which are easily accessible no matter how the case is positioned


#### **SHOCK ABSORPTION**

The iPad is surrounded by rubber bumpers. In addition to the reinforced plastic enclosure this contributes to a carefree protection.

#### **IPAD SUPPORT**

The Rehadapter supports the most recent iPad 10.2", Air, Pro 11" and all upcoming iPads within the same size category and connects hassle-free through the built-in Bluetooth module.

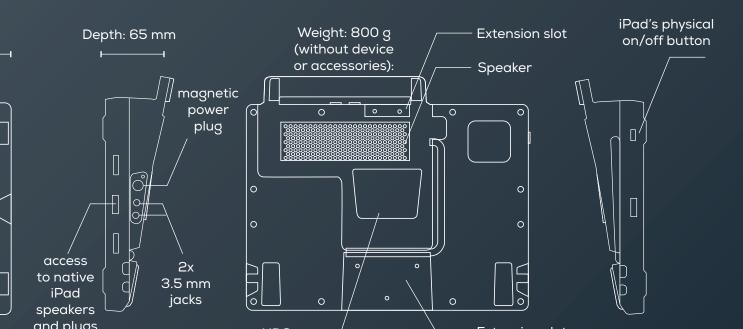

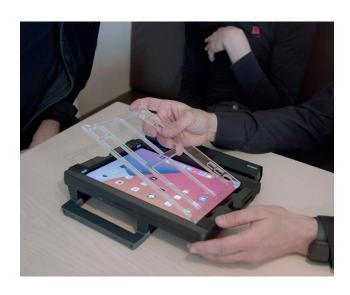



#### **OUR SALES TEAM LOOKS FORWARD TO HEARING FROM YOU**



ROBERT KOHLS
Director of Sales
& Training Europe
rks@rehadapt.com



MAGGIE MAHONEY
CEO Rehadapt
North America
mmy@rehadapt.com







#### **REHADAPT HEADQUARTERS**

Rehadapt Engineering GmbH & Co. KG Heinrich-Hertz-Str. 104 34123 Kassel, Germany

**4** +49 561 220 717 0

**+49 561 220 717 99** 

☑ officeteam@rehadapt.com

#### **REHADAPT NORTH AMERICA**

Rehadapt North America LLC 40 Oser Ave, Suite 40-01 Hauppauge, NY 11788, USA

**4** +1 888 635 6962

☑ officeteam.na@rehadapt.com